


# Quickpin

## Luer Lock Spike

Quickpin is a spike that facilitates easy puncture of small diameter thick rubber stoppers. It provides fast fluid addition and extraction due to the large surface area of the 0.2 micron hydrophobic air-filter that quickly equalizes pressures. It may be used for safe handling of antineoplastic medications.

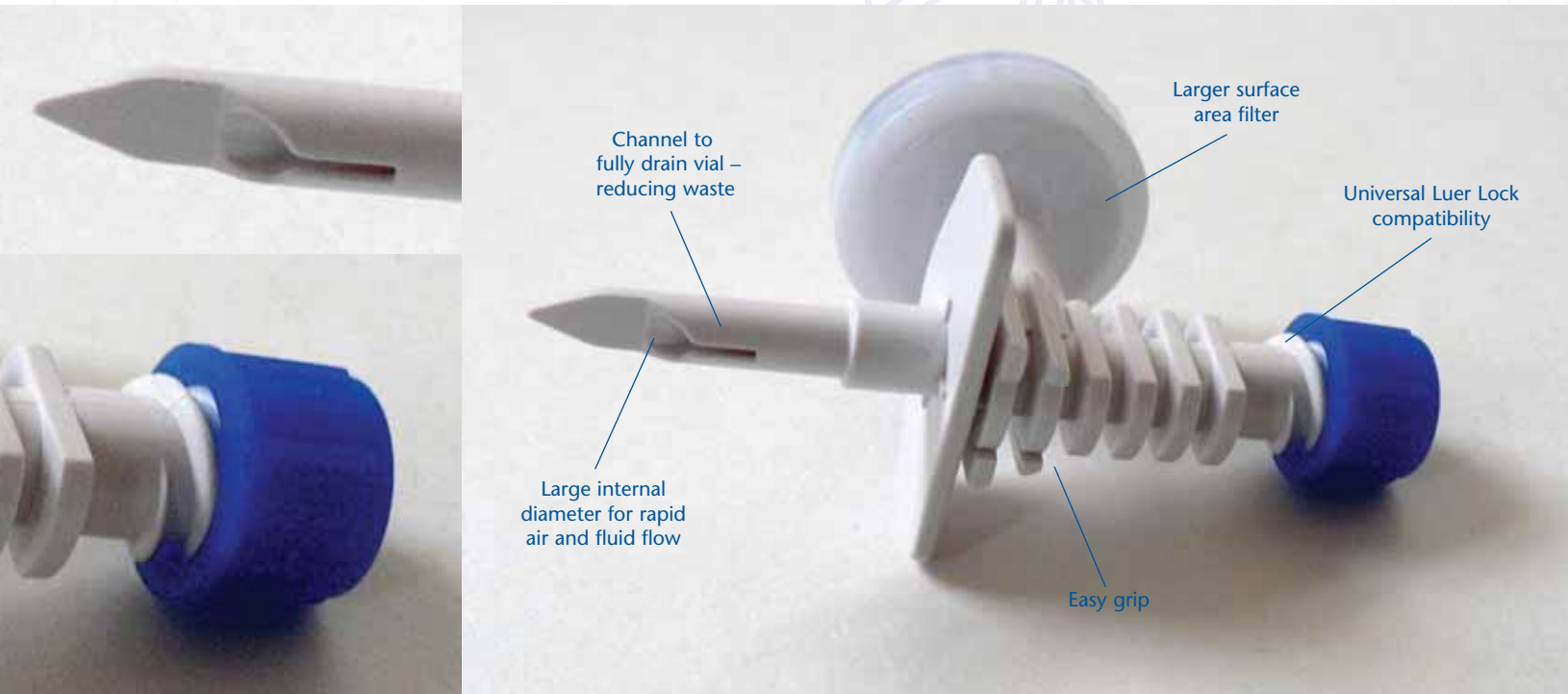

Distributed in USA by:

Grifols USA, LLC

8784 Northwest 18th Terrace  
Doral, FL 33172  
Tel.: [800] 379-0957  
[www.grifolsusa.com](http://www.grifolsusa.com)

- Compatible with oncology agents
- Universal Luer Lock compatibility
- Large internal diameter for rapid air and fluid flow
- Easy grip
- Larger surface area filter that quickly equalizes pressure
- Channel to fully drain vial – reducing waste
- Small priming volume of 0.25mL

Process Control. Protection. Peace of Mind.  
*designed by Grifols*

HD0006-A